


VISCO JET®
- the original

VISCO JET® VJ900

Laboratory Stirrer Set



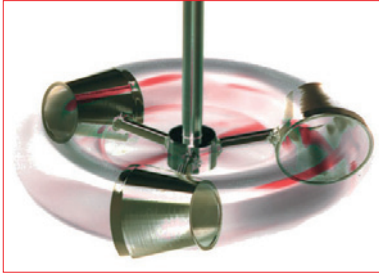

► A quantum leap in mixing technology

The VISCO JET® system is the original one by VISCO JET®, a conical, slow-running, momentum-based range of stirrers. This principle enables the efficient mixing of a range of products even at low circumferential speeds.



► Everything to hand!

The VJ900 Laboratory Stirrer Set contains VISCO JET® stirrer heads in sizes from $d_2=60$ mm to $d_2=150$ mm neatly arranged and immediately to hand for stirring trials in your laboratory. These stirrer heads can be used in all common sizes of vessels from 80 mm to 250 mm in diameter. The shaft contained in the set can be mounted on any laboratory stirrer with a 10 mm chuck.

The spiral version is especially suitable for highly viscous liquids. The VISCO JET® Crack will be used for dispersing of low- and high-viscosities and for stirring in powder materials.



► Description of contents

Case: Silver gray, with aluminum frame extrusion foam insert

Contents: Impeller with M8-thread, patented:

Material 1.4571 (V4A), additional electrolytic polished:

| | | | |
|-------------------|--------|-----------------------|----------------------------------|
| VISCO JET® | 2-fold | $d_2=60$ mm | for vessels Ø approx 80-150 mm* |
| VISCO JET® | 3-fold | $d_2=80$ mm | for vessels Ø approx 120-200 mm* |
| VISCO JET® spiral | 3-fold | $d_2=80$ mm patented | for vessels Ø approx 120-200 mm* |
| VISCO JET® spiral | 3-fold | $d_2=120$ mm patented | for vessels Ø approx 160-250 mm* |
| VISCO JET® CRACK | 3-fold | $d_2=120$ mm | for vessels Ø approx 160-250 mm* |

*depending on the viscosity

Material special plastic POM, M8-brass thread

| | | | |
|---------------|--------|---------------|----------------------------------|
| PROJET® color | 3-fold | $d_2=120$ mm, | for vessels Ø approx 160-250 mm* |
|---------------|--------|---------------|----------------------------------|

*depending on the viscosity

Shaft:

- Diameter 10 mm,
- length 500 mm,
- material stainless steel V4A 1.4571

► Special advantages

- cone principle for gentle stirring
- low rotational speed
- low power requirement
- best stirring results
- no air inclusion
- no foam formation
- no product warming
- also suitable for highly viscous media
- easy to clean

VISCO JET Rührsysteme GmbH

Daimlerstraße 1
D-79761 Waldshut-Tiengen

Phone: +49 7741 96567 0
Fax: +49 7741 96567 15

www.viscojet.com
info@viscojet.com